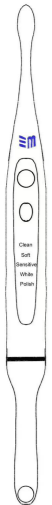
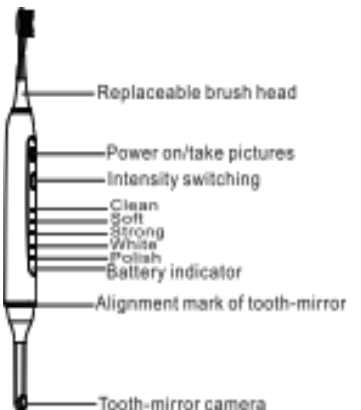
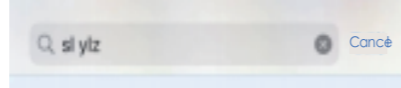
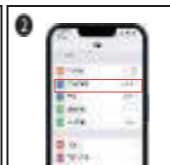
60.00 mm		90.00 mm		<p>Visualized Smart Electric Toothbrush</p> 	<p>1.ProductIntroduction Please be sure to read this introduction manual carefully before using the product in order to use it correctly. This product is a visual smart electric toothbrush with camera. The product features a toothbrush that can be connected to cell phones, tablets and other devices, and the healthy shape of teeth inside your mouth can be seen through APP, and you can take pictures or videos, which provide clear graphical guidance and treatment basis for brushing teeth to remove tooth stains and prevent tooth decay, achieving the effect of whitening teeth and preventing tooth decay from worsening</p> 	<p>2.Safety Precautions To avoid the risk of failure, electric shock, injury, fire and damage to equipment and property, please observe the following safety precautions: (1) This product has a built-in battery. Do not put this product into a fire or expose it to the sun, otherwise the battery may be damaged. (2) Do not use this product when charging. (3) Do not use this product when USB plug is not in place. Only when the plug is properly plugged can it be effectively waterproof. (4) Do not place the tooth-mirror in a humid environment or in liquid to avoid damage caused by immersion in water. Be sure to keep the plug of the tooth-mirror dry. (5) When using the tooth-mirror, make sure that the socket and plug are dry before inserting it, otherwise there will be a risk of damage to the host due to water. (6) Do not modify or repair this product by yourself. (7) When using APP to take pictures as a basis for diagnosis and treatment, be sure to choose a regular dental institution . (8) Do not drop or hit the product hard. (9) Do not charge this product unattended for a long time, the maximum charging time is not more than 12 hours. If you do not use it for a long time, please charge it at least once every six months . (10) Please do not disassemble the product or replace the internal battery by yourself. If you find that the internal battery has swelling, leakage or other abnormalities, please stop using it immediately and contact our company's after-service .</p>	<p>TM04 visual APP installation and use tutorial 1. Please search for sl ylz in App Store/Google play to download and install.</p>    <p>Scan the QR code with your phone, download the APP(slylz) and install it.</p>	<p>Android China Android Google Play IOS app store</p> <p>IOS After scanning the code, open the IOS App Store to download and install. You can also search for SL ylz to download and install in the App Store directly</p> <p>Android After scanning the code, open the Android google play to download and install it through a third-party browser.</p> <p>Android China Android google play IOS app store</p> <p>Attention: If there is a safety reminder during the installation process, choose to continue the installation</p>	<p>2. APP connection instructions</p>       <p>Return to the app, "Connected" is displayed, Click "Connected."</p> <p>Enter the camera screen, as shown above,</p>	<p>*Charging instructions 1. Press and hold the button to power off and then charge 2. The red light blinks when charging. 3. When fully charged, the green light keeps on .</p> <p>* Troubleshooting 1. The picture is stuck A. The data transmission distance is too far, do not stay away from the phone when using the product. B. Insufficient power will also lead to stalling, please charge in time . 2. WiFi signal cannot be transmitted/intermittently A. The lower the power, the weaker the WiFi signal, fully charged and then used. B. Over the WiFi signal transmission distance, the closer the product is to the phone, the better. C. There are too many WiFi signals around, it will affect the signal of the product, try to use in the area with less WiFi signal 3. No image is displayed A. Check whether the phone is connected to the correct product WIFI name . B. After confirming the connection, if there is no display, close the APP and product, open it again and connect again .</p>	<p>* Teeth Brushing Step one Dampen the bristles and apply some toothpaste to the bristles.</p>  <p>Step two Press the on/off button to turn on the toothbrush. Power on/Off/Mode switch In special cases, you can press and hold for more than 5 seconds to power off. The brushing process should be gently applied a little pressure to maximize the effect of the toothbrush. Slowly move the toothbrush at a 45° Angle between the bristles and gums. During the brushing process, press the button to switch the brushing mode freely. The recommended brushing time is 2 minutes.</p> 	<p>* Clean brush head After brushing your teeth, clean the toothbrush head with water (Be sure to cover the charging bottom when cleaning).</p>  <p>* Maintenance instructions 1. Wash the brush head and bristles after each use, and air dry naturally (do not pull the bristles, which is easy to deform and fall off). 2. Wipe with a dry cloth after each use and air dry naturally. 3. Do not use water higher than 40C or acid cleaning solution to clean the brush handle, charging cable, power cord and adapter, the adapter is strictly prohibited in water. 4. Before cleaning the charging device, be sure to pull out the adapter from the socket, otherwise it may cause electric shock. 5. It is recommended to replace the brush head every three months, and replace it when the bristles are bent or deformed.</p> <p>* Environmental protection - When disposing of products, please do not place them with the general domestic garbage. - Built-in rechargeable batteries contain substances that may pollute the environment. When discarding products, do not throw them away randomly. Please place the products in front of the government designated recycling station.</p>	<p>* Warranty card Thank you for your support and trust ! Electric toothbrush should be able to play its advantages under normal use and maintenance, and any damage caused by manufacturing process or components within one year of purchase will be replaced free of charge.</p> <p>* Warranty Terms 1. Free warranty period, starting from the date of purchase. 2. Warranty products, please show the following voucher or purchase record. 3. Replacement of spare parts more than one year, according to the replacement parts price charges. 4. Toothbrush head is easy to wear products, not within the scope of warranty. 5. Under any of the following circumstances, free warranty is not implemented. 1) Remove the fuselage during use, causing the toothbrush to not be used normally . 2) Unable to provide purchase vouchers or purchase records</p>	<p>Product model _____ Date of Purchase _____ Place of Purchase _____ Customer Name _____ Mobile phone _____ Contact Address _____ Fault Description _____</p> <p>* Argument Model: TM04 Rated input: DC5V 1A Operating voltage: 3.0V-4.2V Battery capacity: 1500 mah Power: 2w Waterproof class: IPX7 Company: Shenzhen THREEMAN Technology Co., LTD Address: Building 4, YingJinNianHua Garden, NanWan Street, NanLing, LongGang District, Shenzhen</p> 
----------	--	----------	--	---	--	--	---	--	--	--	--	--	---	--

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